

LDD 100



Defined dilution system for large droplets up to 10 m

Model Variations



LDD 100 H

Heatable dilution system up to 150 °C for large droplets up to 10 m

Description

The dilution of large droplets is particularly important when measuring highly concentrated droplet aerosols. Since large droplets are difficult to dilute, standard systems only work up to a size of 1 - 2 m. The newly developed dilution system LDD 100 (dilution factor 100) is the first system to dilute almost loss-free large droplets up to 10 m.

The good dilution factor of large droplets was tested with monodisperse DEHS droplets (oil) of different sizes. The results for the sizes 5 m and 7 m are shown in Table 1.

LDD 100:

Particle size	Number counts without dilution	Number counts with dilution	Dilution factor
5 μm	304322	3043	100.01
7 μm	236687	2370	99.87

Chart 1: Dilution of monodisperse DEHS droplets with LDD 100

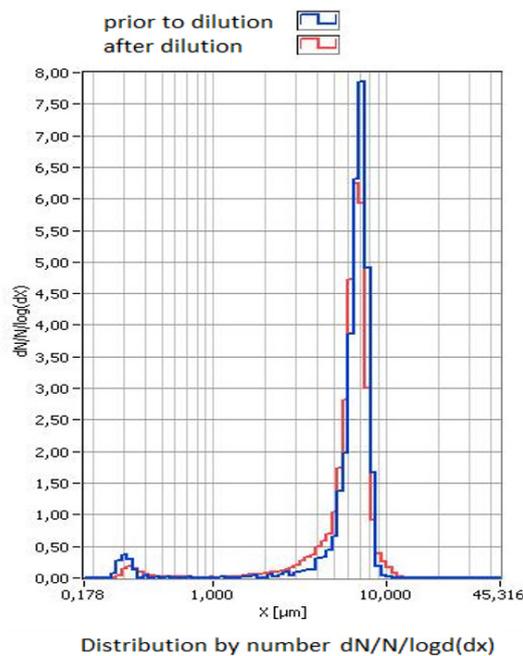


Fig. 1: Distribution of LDD 100 (7 m)

Benefits

- Defined dilution of large droplets of factor 100
- Proven dilution factor 100 for droplet sizes up to 7 μ m
- Easy connection with Promo[®] and welas[®] digital aerosol spectrometers
- Internal pump for autonomous operations
- Resistant to pressure fluctuations of \pm 200 mbar
- Simple handling
- Robust, durable, low maintenance
- Cost effective

Applications

- Measurement of blow-by aerosols according to ISO 17536
- Dilution of compressed air
- Measurement of cooling lubricant aerosols

Palas GmbH
Partikel- und Lasermesstechnik
Greschbachstrasse 3 b
76229 Karlsruhe
Germany

Managing Partner:
Dr.-Ing. Maximilian Weiß, Udo Fuchslocher
Commercial Register:
register court: Mannheim
company registration number: HRB 103813
USt-Id: DE143585902



Contact: E-Mail: mail@palas.de Internet: www.palas.de Tel: +49 (0)721 96213-0 Fax: +49 (0)721 96213-33